concept

**Sonický zubní kartáček
& Přístroj na mezizubní hygienu**
**Elektrická sonická zubná kefka
& Prístroj na medzizubnú hygienu**
Szczoteczka soniczna & Irygator
**Szónikus elektromos fogkefe
& Fogköztisztító**
**Sonic elektriskā zobu birste
& Zobstarpu tīrītājs**
**Sonic electric toothbrush
& Interdental Cleaner**

**Elektrische Schallzahnbürste
& Gerät zur Zahnzwischenraumreinigung**
**Brosse à dents électrique
& Microjet interdentaire**
**Spazzolino sonico ricaricabile
& Sistema di pulizia interdentale**
**Cepillo dental eléctrico sónico
& Limpiador interdental**
**Periuță de dinți sonică electrică
& Aparat de curățare interdentară**

ZK4030

CZ

SK

PL

HU

LV

EN

DE

FR

IT

ES

RO

PODĚKOVÁNÍ

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ACKNOWLEDGEMENTS

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

Technical specifications	
Voltage	100-240 V ~ 50/60 Hz 5 V / DC
Input	1.8 W
Function	5 cleaning modes (Clean/Sensitive/White/Gum Care/Massage) 3 operating modes (Normal/Soft/Pulse)
Charging	wireless induction charger
Water resistance	IPX7

IMPORTANT SAFETY PRECAUTIONS

- Do not use the device in any other way than as described in this manual.
- If you have had dental surgery in the last two months, consult your dentist about using the device.
- Do not use the device if you have a mouth injury.
- Do not use attachments other than those recommended by the manufacturer.
- The device is intended for dental hygiene only. Do not use it for other purposes.
- If you experience discomfort or pain, stop using the device and seek medical attention.
- Remove all the covering and marketing materials from the device before the first use.
- Make sure the mains voltage corresponds to the values on the rating plate of the device.
- Do not leave the device unattended when it is switched on or connected to the mains.
- If you put the device down during use, always switch it off.
- Before connecting or disconnecting the device from the mains socket, make sure that the device is switched off.
- When disconnecting the device from the wall outlet, never pull the power cord, grasp the plug end and disconnect it by removing it from the receptacle.

- Do not allow children or unauthorized people to handle the device, use it out of their reach.
- Do not use the device in the bathtub.
- Do not leave the device in a place from which it can fall into a sink or a bathtub.
- If the device has fallen into water, do not reach for it and unplug it immediately.
- Do not share the attachments with other people.
- Do not place metal objects (coins, staples, etc.) on the charger.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the device only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the device is used near children.
- Do not use the device as a toy.
- **Before cleaning and after use of the device, switch the device off and disconnect it from the mains.**
- Do not let the power cord touch hot surfaces.
- Do not allow the device to hang on the power cord.
- Do not use the device in an environment where there are explosive gases and flammable substances (solvents, varnishes, adhesives etc.).
- Do not use coarse and chemically aggressive substances to clean the device.
- Do not use the device if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- Do not use the device with a damaged cord or plug and have it repaired immediately by an authorized service centre.
- Do not wrap the power cord around the device body.
- Do not use the device outdoors.
- The device is only suitable for home use and is not intended for commercial use.
- **Do not immerse the power cord, plug or device in water or any other liquid.**
- Do not use any accessories other than those recommended by the manufacturer.
- Do not repair the device yourself. Contact an authorized service centre.

- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the device in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are 8 years old and under supervision. Children under the age of 8 must be kept away from the device and its power supply. Do not let children play with the device.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.


CAUTION:

This device is waterproof. It is suitable for cleaning under running water. Do not use or clean the device while it is connected to the power cord.

PRODUCT DESCRIPTION

- 1 Power button / Mode switch button
- 2 Clean, Sensitive, White, Gum Care, Massage modes
- 3 Brush for regular cleaning
- 4 Power button
- 5 Normal, Soft, Pulse modes
- 6 Brush with fine bristles
- 7 Attachments - nozzles
- 8 Wireless induction charger

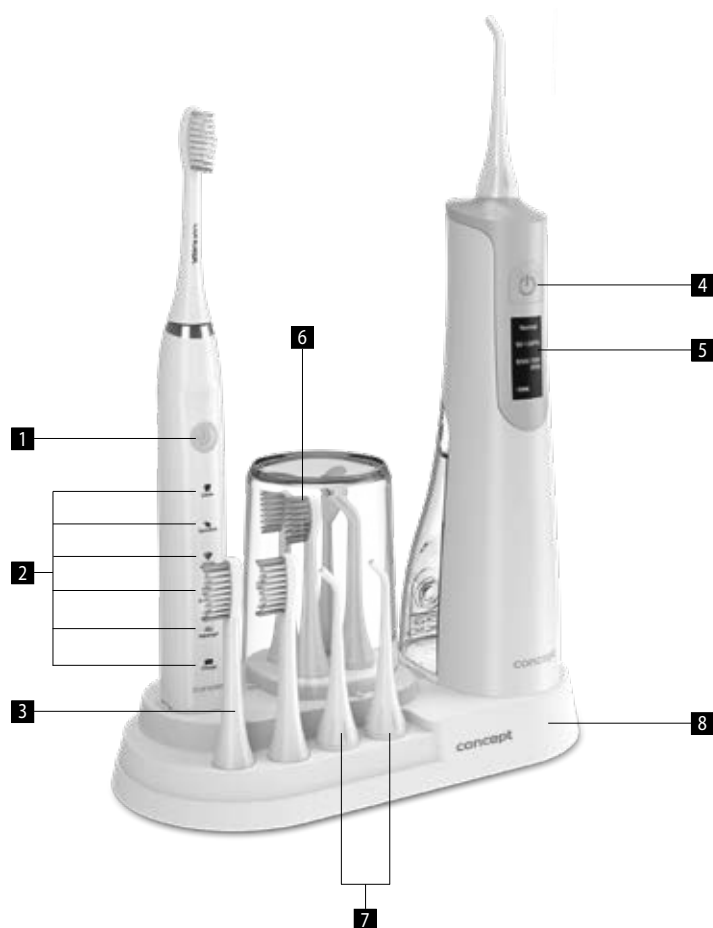OPERATING INSTRUCTIONS

Charging

- Fully charge the battery before using the device for the first time.
- Charge only when the battery is low.
- To make better use of the battery, we recommend discharging the battery once a year.

SONIC ELECTRIC TOOTHBRUSH

Modes

Clean - standard choice for cleaning, removes 100 % of plaque

Sensitive - gentle cleaning of sensitive teeth and gums


White - helps to remove surface stains

Gum Care - gum cleaning

Massage - pulsating vibrations for a gentle gum massage

Smart timer


- The recommended cleaning time is at least 2 minutes. The device switches off automatically after two minutes.
- Every 30 seconds, the device pauses to alert you to move to the next part of your mouth.

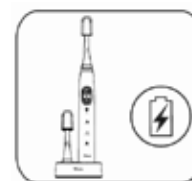Fully charge the battery before using the device for the first time - the indicator lights up green.

To clean your gums effectively, perform the cleaning at an angle of 45°.

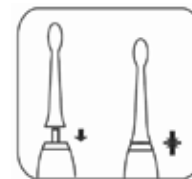
Place the brush head on the brush body.

Apply the brush lightly to avoid damaging the gums.

Rinse the brush head, apply the paste, put the brush head in your mouth and then turn on the brush.

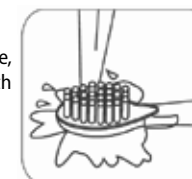Rinse the brush head with water and dry.

INTERDENTAL CLEANER

Preparation

- Fill the tank with water and close it tightly.
- Use tap water or heated water.
- Do not use salt water, toothpaste or chemicals.
- Attach the nozzle.

If the device is switched on but the water does not spray out:

- Make sure there are no foreign objects in the nozzle.
- Fill the tank with water.
- Turn on the device and carefully turn it upside down several times OR remove the nozzle, turn on the device and pour water into the nozzle socket.

Use

- Insert the nozzle into your mouth towards your teeth and lean over the sink.
- Turn on the device.
- Hold the device vertically with the nozzle pointing up.
- Close your mouth so that the water does not spray out.
- Move smoothly along the teeth.
- A full tank lasts for approx. 45 seconds of use.

Modes

Normal - regular cleaning

Soft - gentle cleaning of teeth and gums

Pulse - effective removal of residues and plaque

To set the mode, press the power button. To stop the water flow, press the button for 5 seconds.

Nozzle types

Basic cleaning nozzle - removes residues and plaque

Nozzle for teeth with braces - cleans and rinses hard-to-reach areas around braces

Nozzle for cleaning the tongue - removes odours and keeps breath fresh, use "pulse" mode and gently run the nozzle over the tongue

TROUBLE SHOOTING**Sonic electric toothbrush**

Problem	Possible cause	Remedy
Weak vibrations	Too much pressure on teeth.	Reduce pressure on teeth.
	Improperly fitted brush head	Fit the brush head again
Loose brush head or abnormal sound	The brush head is damaged	Replace the brush head
Shorter operating time with a charged battery	Charging error.	Charge the toothbrush correctly.
	Empty battery	Charge the brush to full charge
The toothbrush does not work	Power has run out	Charge the battery

Interdental cleaner

Problem	Possible cause	Remedy
The device only works for a few minutes (even after charging)	End of battery life.	Contact customer service
	Insufficient charging time	Charge the device in one go for more than 16 hours
Water pressure too low	The nozzle is dirty/deformed	Replace the nozzle
The water does not spray out	The device is too tilted during use	see Preparation above
	There is dirt in the nozzle	
	The nozzle is incorrectly fitted	
The device does not work	Newly purchased product or unused product for more than three months	Charge the device for at least 10 hours

CLEANING AND MAINTENANCE**Sonic electric toothbrush**

Brush head - After use, rinse with water and dry.

Brush body - Rinse with water and dry. Do not use water warmer than 40 °C or any cleaning agents.

Charging base - Wipe only with a dry cloth. Wet and damp fabrics can cause a short circuit.

Replacing the brush head - Dentists recommend replacing the brush head every 3 months.

Store the device in a dry and dust-free place.

Do not expose to high temperatures.

Keep out of reach of children.

Interdental cleaner

1. Pour water out of the tank.
 - Switch on the device and let the remaining water drain out.
 - If you have used a liquid other than water, wash the tank thoroughly. Then fill the tank with warm water, point the nozzle at the sink and allow all the water to drain. This prevents clogging of the nozzle.
2. Press the nozzle release button, remove the nozzle and place it in the case.
3. Clean the device with a soft cloth or a mild, non-abrasive cleaning agent.

Nozzles

- Press the nozzle release button and remove the nozzle.
- Rinse with water and wipe dry.
- We recommend replacing the nozzle after 6 months. If the nozzle is deformed, replace it immediately.

Device body

- Wipe with a cloth and a mild detergent. Wipe dry.
- In case of heavy soiling, rinse with water. Then dry.
- Do not use water warmer than 50 °C.

Tank

- Rinse with water and wipe with a cloth.
- If you do not use the device for a prolonged period of time, let the tank dry.

Use a mild detergent to clean. DO NOT use gasoline or alcohol.

SERVICE CENTRE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified specialist or the service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan:

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (waste electrical and electronic equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without prior notice and we reserve the right to change them.

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Warranty

The manufacturer (or importer) is responsible for ensuring the product complies with the requirements of applicable legal regulations as well as those of the relevant technical standards. Moreover, they are responsible for ensuring the product has the properties the manufacturer described in documents related to the goods or those reasonably expected by the customer with regard to the nature of the goods or based on advertising produced by the manufacturer, and further they are responsible for ensuring the product is fit for the purpose proposed by the manufacturer or that a product of the same type is normally used for.

The quality warranty term is 24 months from product takeover by the customer.

The warranty does not apply to wear and tear caused by regular use. The customer shall not be entitled to any warranty claims if, prior to taking the product over, they knew the product contained a defect or if the defect is attributable to the customer.

In terms of these warranty conditions, a defective battery is considered a battery that is electronically inoperative (i.e., does not communicate with the base unit, the charging or discharging circuits are defective, etc.) or has an internal short circuit.

Defects or changes in parameters caused by normal use and wear (e.g. natural capacity decrease) are in no way considered a defect in terms of warranty.

The battery is provided warranty service for a period of 6 months. In terms of these warranty conditions, the battery has the nature of a consumable, for more details see „Consumables“.

Consumables

For consumables (batteries, etc.) we guarantee the preservation of technical parameters for 6 months, in case it was not damaged by improper use and its lifetime has not expired.

The warranty specifically does not apply:

- if the product installation, operation and service conditions stipulated in the product operating manual have not been adhered to,
- to malfunctions caused due to mechanical, heat or chemical damage, short circuit, over voltage or incorrect installation,
- to malfunctions caused by an inexpert third-party intervention,
- to malfunctions caused by natural disaster,
- to malfunctions caused by insufficient or inappropriate maintenance in violation of the operating manual, including malfunctions caused by water and other sediments,
- to changes in colour of the heating surface or to scratching of the surface caused as a result of using the products in an unusual manner,
- to appearance and functional changes caused by exposure to sunlight, thermal radiation of water and other sediments,
- if the service life of certain product parts expires, e.g. for accumulators, bulbs, etc.

The warranty does not apply to any products and services provided along with the product (gifts, promotional articles, etc.).

Filing a complaint

A complaint against a product defect must be filed as soon as identified, yet no later than before the end of the warranty term.

The customer must file a product complaint with the dealer from which they have purchased the product, or with any authorised service centre, a list of which is included in the product package, or available at www.my-concept.com.

While filing a product complaint, the product must be duly cleaned and securely packed so as to prevent any damage during its transport to an authorised service centre, where relevant, unless the product is delivered in person.

The customer must submit proof of having concluded a purchase contract for the product by producing the receipt.

While filing their complaint, the customer must indicate the noted defect and identify the preferred complaint application method.

Complaints processing

As long as the noted defect may be removed, the user has the right to have the defect duly removed free of charge on a timely basis.

If the noted defect is not removable, or if the customer becomes entitled to replacement of the product or a part thereof, yet the replacement is not possible, for example due to the product having been sold out, the customer has the right to return the product (withdrawal from the contract).

The customer shall have the right to a new product (replacement) or to replacement of a part of the product even if the defect can be removed, provided they cannot properly use the product due to repeated occurrence of the defect or due to a high number of such defects. In such a case, the customer also has the right to return the product (by withdrawing from the contract).

If the product is not returned (the customer does not withdraw from the contract), or if the customer does not apply the right to a new defect-free product (replacement), to replacement of a part thereof or to repair of the product, they may request a reasonable discount. The customer also has a right to a reasonable discount if a new defect-free product cannot be supplied to them, or if a product part cannot be replaced or the product repaired unless the situation is remedied within a reasonable time limit, or if remedying the situation would create major discomfort on the part of the customer.

The seller, authorised service centre or a staff member authorised by them must decide about each complaint immediately or within three business days in complicated cases. This term does not include a reasonable period of time, depending on the type of product concerned, required for the defect to be assessed by an expert. A complaint, including defect removal, must be dealt with without any undue delay, yet no later than within 30 calendar days of the complaint filing date, unless the seller and the customer agree on a later deadline.

On returning the product (withdrawing from the contract) the customer must return any accessories and documents supplied along with the product.

The customer does not have the right to keep the defective parts and components of the product replaced as part of a repair of the product.

This shall be without prejudice to any other rights the customer may have in relation to the purchase of the product.

Remark: Complaints against products damaged in transport are governed by the carrier's complaints procedure.

Manufacturer

Jindřich Valenta - Concept
Vysokomyšská 1800, 565 01 Choceň
Czech Republic
tel.: +420 465 471 400
fax +420 465 473 304
Company ID No. 13216660
email: servis@my-concept.cz
www: [www: www.my-concept.com](http://www.my-concept.com)

Product data

Model:	
Production number:	
Date of purchase:	Seal and signature of vendor:

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Verners VT Ltd.