# niceboy ion®

## SONIC **KIDS**

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
**PACKAGE CONTAINS**


Sonic Niceboy Kids toothbrush handle  
 1x brush heads (soft, medium, hard)  
 Wireless docking station  
 Micro charging cable USB

1x brush heads

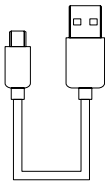Charging USB cable

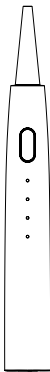body of sonic toothbrush

**INSTRUCTIONS**

1. Fully charge the toothbrush before first use. The charging connector can be accessed by removing the cap on the bottom of the toothbrush.


**The charging process:**

- 0 % -10 %, charging indicator shows red
- 10-100 %, charging indicator shows orange
- Once charging is complete, the charging indicator shows green


*Note:*

If the product has not been used for a long time, the charging indicator may not be displayed when charging the electric toothbrush. To see the charging status, wait patiently for 3-5 minutes.


2. Insert the new brush head as shown; if you need to remove it, just pull it out as shown.

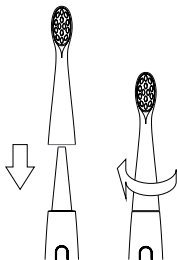3. Wet the bristles and squeeze out a reasonable amount of toothpaste. Touch the toothpaste-coated bristles to the surface of the tooth in any part of the mouth.
4. Press the switch and move the jaw up and down slightly to clean your teeth from the outside and inside. The oral cavity is divided into 4 areas and the recommended cleaning time for each area is 30 seconds. When you have finished cleaning that area for 30 seconds, the electric toothbrush will stop briefly to remind you to clean the next area.

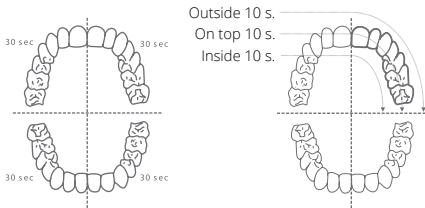
5. Press the mode switch button on the handle to switch the mode according to the person's preference.

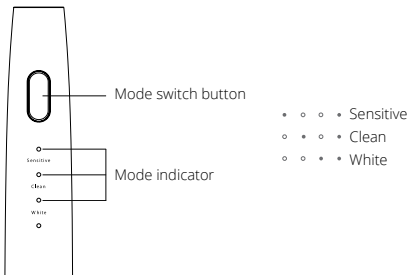*Note:*

1. Sensitive mode is activated by default
2. Switching mode: within 5 seconds after switching on, briefly press the mode switch button to switch modes
3. Auto memory mode: use any mode for more than 5 seconds, the system will record this mode. The memory mode is activated by default the next time the power is turned on.
6. After the set toothbrushing time of 2 minutes, the electric toothbrush automatically stops. After brushing your teeth, rinse the brush head with running water and shake off any excess water from the brush head and handle. Keep the toothbrush crystals dry and store in a dry place.

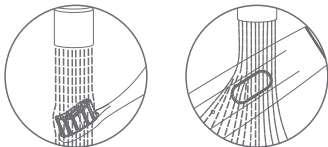*Note:*

1. If you want to stop the electric toothbrush in the middle of cleaning your mouth, press the switch briefly for 5 seconds after switching on the electric toothbrush.
2. If you are not going to use it for a long time, disconnect the charging plug and clean the brush head and handle thoroughly, store in a cool and dry place
3. It is recommended to replace the toothbrush head every 60 days.

## **SAFETY INSTRUCTIONS**

Read this manual carefully before using the device and keep it for future reference. Keep the connected charging cable out of the water. Do not place or store it in or near a filled bathtub, sink, etc. Do not immerse the cable in water or any other liquid. After cleaning, make sure the cable is completely dry before connecting it to the power supply. Before connecting the appliance, check that the voltage indicated on the appliance corresponds to that of the local power supply. Never use the charging cable if it is damaged in any way. If the cable or connector is damaged, always have it replaced with a new one to avoid danger. If the device is damaged in any way (brush head, toothbrush handle), stop using it. This device does not contain any parts that can be repaired by the user. If the appliance is damaged, contact the complaints centre. Do not charge the appliance near heated surfaces. This product is not intended for use by persons (including children) with reduced physical, sensory or mental abilities or lack of experience and knowledge, unless they have been supervised or instructed in the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. Improper handling (prolonged charging, short-circuiting, breakage by another object, etc.) can lead to fires, overheating or leakage of the battery, for example. The battery or accumulator that may be included in the product is subject to a lifetime of six months as it is a consumable material.

## **PRECAUTIONARY MEASURES**

Do not use the toothbrush in the bath or shower. If you have had gum or oral surgery in the previous 2 months, consult your dentist before using the toothbrush. If excessive bleeding occurs after using this toothbrush or if bleeding continues after 1 week of use, consult your dentist. The Niceboy ION Sonic Kids toothbrush meets the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, consult your doctor or the manufacturer of the implanted device before using the Niceboy ION Sonic Kids toothbrush. This device was designed for cleaning teeth, gums and tongue only. Do not use it for any other purpose. If you experience any discomfort or pain, stop using the device and contact your doctor. If you have a medical condition, consult your physician before using the Niceboy ION Sonic Kids. The Niceboy ION Sonic Kids toothbrush is a personal care device and is not intended for use with multiple patients in a dental office. Do not use a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not clean the brush head in the dishwasher. Do not use brush heads other than those recommended by the manufacturer. Supervise children when brushing their teeth until they have developed good brushing habits. If the toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), clean the brush head thoroughly with soap and water after each use. This will prevent possible cracking of the plastic.

## PRODUCT SPECIFICATIONS

Rotational speed:	60 000
Voltage / Power:	3,7V == 1,2W
Charging:	Standard Micro USB charger
Output Voltage:	5V == 300 mA
Operating time:	40 days (calculated based on operating condition of 2 minutes/time and 2 times/day)
Battery capacity:	500 mAh
Charging time:	2,5 hours
Waterproof:	IPX7

NICEBOY s.r.o. hereby declares that the type of radio equipment NICEBOY ION SONIC KIDS complies with Directives 2014/53 / EU, 2014/30 / EU, 2014/35 / EU, and 2011/65 / EU. The full content of EU Declaration of Conformity is available on the following websites: <https://niceboy.eu/sk/declaration/ion-sonic-kids>

## INFORMATION FOR USERS ON THE DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT (HOUSEHOLDS)

The symbol on the product or in the accompanying documentation means that used electrical or electronic products must not be disposed of with municipal waste. In order to dispose of the product correctly, please take it to the designated collection points where it will be accepted free of charge. Proper disposal of this product will help conserve valuable natural resources and help prevent potential negative environmental and human health impacts that could result from improper waste disposal. For more details, contact your local authority or nearest collection point. Fines may be imposed for improper disposal of this type of waste in accordance with national regulations.

Information for users on the disposal of electrical and electronic devices (corporate and business use)

For proper disposal of electrical and electronic equipment for corporate and business use, consult the manufacturer or importer of this product. They will provide you with information on how to dispose of the product and, depending on the date of placing the electrical equipment on the market, who is obliged to finance the disposal of this electrical equipment. Information on disposal in other countries outside the European Union. The above symbol is only valid in European Union countries. For the correct disposal of electrical and electronic equipment, please ask your authorities or equipment dealer for detailed information.

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